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
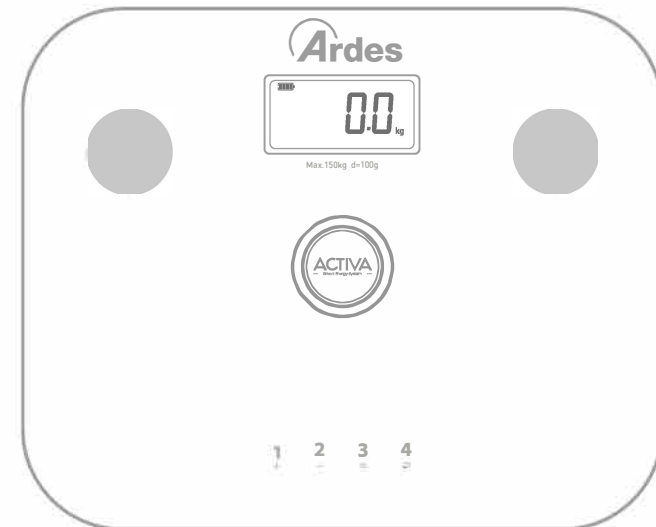
**Thank you for your purchase.**

For the proper use and maintenance of the scale, please read these instructions completely before operating.


## USER GUIDE
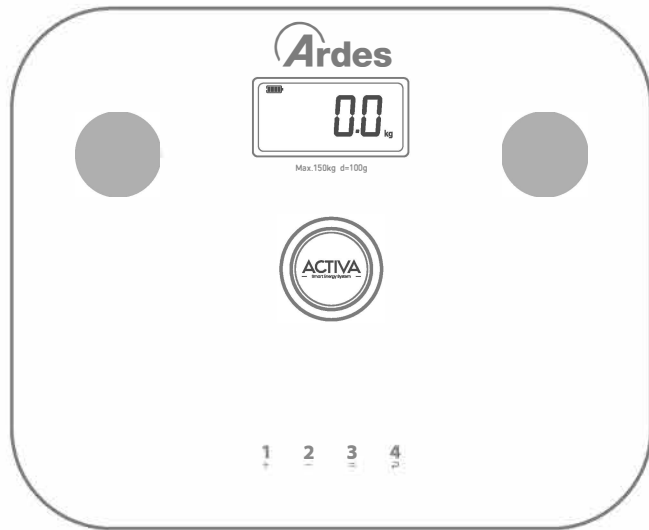
## BODY INGREDIENT ANALYSIS SCALE

AR2PPBFM

### Key function

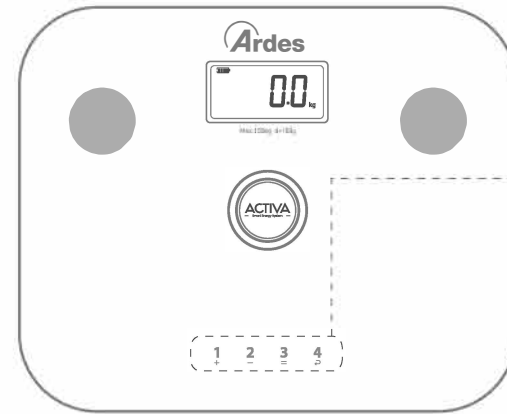• Touch Keys

**1** +      **2** -      **3** =      **4** ↺

Press **1** to increase the data or switch the gender.


Press **2** to reduce the data or switch the gender.

Press **3** to confirm the data input.

Press **4** to go back to the last step. Hold down the button if you want to adjust the settings.

### How to use

Before using, put the scale on a hard and flat surface. There are two modes to start.
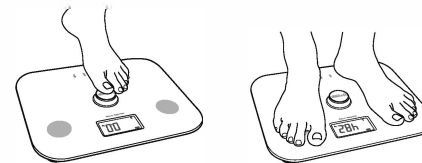
#### Mode 1: For weighing only, no ingredient analysis

- At least two quick presses on the ACTIVA button will generate power for weighing. Please wait until "0.0" appears.

⚠ Do not step on the scale until "0.0" appears.

#### Mode 2: For weighing and ingredient analysis

**Step 1:** At least two quick presses on the power button by your feet will turn it on.

Family user number selection—press and hold down any of the keys **1** / **2** / **3** / **4** to set and input your related personal data into the memory before using for the first time.


#### Step 2: Gender Selection

When the "MALE" is twinkling

Press **1** or **2** to select male or female

Press **3** to confirm
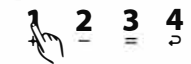
#### Step 3: Units setting :

When the unit is twinkling

Press **1** or **2** to change unit

Press **3** to confirm


**Step 4 : Height setting :**  
When the height number is twinkling  
Press **1** or **2** to change height  
Press **3** to confirm

**Step 4 : Age setting :**  
When the age number is twinkling  
Press **1** or **2** to change age  
Press **3** to confirm  
Press one of the buttons for a long time if you want to adjust the parameter quickly.

**Step 5: Weighing and measuring:**

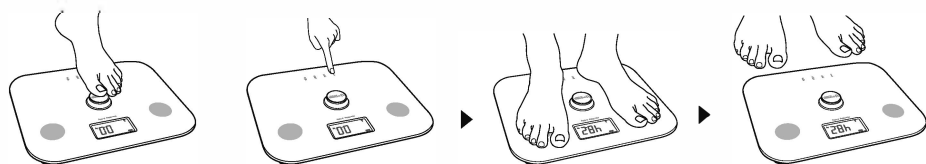When "0..." is displayed from left side to right side, it is analysing.  
Please do not leave the scale until your weight is locked on the screen.

- NOTE**
- 1)When "0.0" appears, stand with clean, dry bare feet on the scale. Your feet must stand evenly on the electrodes. Do not stand outside of the electrodes.
  - 2)When using for the second time, you can begin by pressing **1/2/3/4** for few seconds to begin measuring directly.
  - 3)This button "**UNIT**" can not change the units, please change the units on the Step 2 of mode 2.

## Manual

### Error Indication

**Over-load Indication**  
When the scale is over-loaded[more than 150kg], the screen will appear "Err".

It means the wrong operation caused the fat percentage Err2.

**Low battery indication**  
When the screen indication "**☐**", it will turn off automatically after few seconds. Simply reload it by stamping on the U-power button!

## Specification

Equipped with a high precision strain gauge sensor system  
Max. Weight: 150kgs/330lbs  
Min. Effective : 3kg/6.61b  
Units:KG/LB  
Body fat/Muscle/Water division: 0.1%  
Body bone division: 0.1kg  
BMI: 0.1  
Age range for body fat: 10-80 years  
Height range: 100-220cm/3ft-3in-7ft-3in  
Body bone range: 0-15kg

Muscle range: 13.0-100.0%  
Fat range: 5-50%  
Water range: 35-75%  
Height unit : cm/ft  
Memory: 4 persons  
Over-load / Low power indication  
Auto zero resetting / Auto power off  
Power : Automatic power generation control  
Battery-free technology

## Attention

### Attention for safety

- The scale should not be used for anyone who has a pacemaker, or with any other medical transplanting. The impedance passed through our body will have disturbance of the normal working of medical equipment.
- This scale is for family use only. We do not advise you to use it as medical equipment.
- Do not use it on a wet or slippery floor. Do not stand at the edge of the scale to avoid falling off and injury.
- We advise pregnant women, elderly persons and disabled persons to use the scale with help.

### Attention for operation

**The body fat reading may vary depending on the amount of body water present. The body fat percentage reading may vary for the following reasons: diet, alcohol consumption, heavy exercise, certain prescription drugs or diuretics, illness, or a woman's menstrual cycle.**

- Please take off your socks, and keep bare feet clean and dry when using.
- We advise you use the scale in bare feet. Wearing socks or shoes will impede the accuracy of the reading.
- Extreme water loss and dieting will have a negative effect on the reading.
- Heavy exercise will have a negative effect on the reading.
- An unstable or incorrect reading can be caused by not connecting with the metal sheet fully. Please clean dirt off the scale and step on the scale as per instructions.
- Keep stable when using. Any shaking will lead to an inaccurate reading.

The scale should not be used by the following groups:

1. People under 10 years old, or over 80 years old.
  2. Athletic or muscular body type (body builders).
  3. Pregnant women.
- For the above group, the body scale reading is only for reference.

## Maintenance

**Please place the scale on a clean and level surface.**

01. Please do not use the scale in extreme high or low temperature environment.
02. Please do not jump or stamp on the scale and avoid dropping it, as the scale is a highly precise electronic instrument.
03. To prevent scratches on glass panel, it is advised not to wear any footwear while taking any measurements.
04. When cleaning the scale, please use only normal cloth and no chemical cleaners should be added to avoid any damage to the scale.
05. This scale is only designed for home use and should not be used commercially for any reason at all. The < Body Fat Ratio Chart > is for reference purpose only.
06. Please remove battery from the scale if it is not used frequently.

## Frequently asked questions

### Is there anything wrong with "----" displayed on the screen?

- The scale is auto zero setting, please wait for a second, "0" will be displayed.
- No parameter was found in scale, please input your personal parameter.

### The scale only displayed weight, and not any body fat percentage is found?

When measuring the body fat, shaking will disturb the body fat percentage measuring. Please keep your body straight when weighing, and be sure that your foot is bare and keep in contact with metal sheet on the scale.

### The body fat percentage is much higher than normal?

Please double-check your personal parameter covering gender, height and age. Make sure your bare foot is clean and dry, and is in contact with the scale. Also, any major change or loss in water will affect the reading of body fat scale. Make sure you stand still when taking your reading.

## Schedule for reference

### Body fat percentage ( % )

Age	Female				Male			
	Underfat(-)	Healthy(0)	Overfat(+)	Obese(++)	Underfat(-)	Healthy(0)	Overfat(+)	Obese(++)
10-20	≤ 19.5%	19.6-33.5%	33.6-38.5%	≧ 38.6%	≤ 9.5%	9.6-20.5%	20.6-25.5%	≧ 25.6%
21-30	≤ 20.0%	20.1-34.0%	34.1-39.0%	≧ 39.1%	≤ 10.0%	10.1-21.0%	21.1-26.0%	≧ 26.1%
31-40	≤ 20.5%	20.6-34.5%	34.6-39.5%	≧ 39.6%	≤ 10.5%	10.6-21.5%	21.6-26.5%	≧ 26.6%
41-50	≤ 21.0%	21.1-35.0%	35.1-40.0%	≧ 40.1%	≤ 11.0%	11.1-22.0%	22.1-27.0%	≧ 27.1%
51-60	≤ 21.5%	21.6-35.5%	35.6-40.5%	≧ 40.6%	≤ 12.0%	12.1-23.0%	23.1-28.0%	≧ 28.1%
61-70	≤ 22.0%	22.1-36.0%	36.1-41.0%	≧ 41.1%	≤ 13.0%	13.1-24.0%	24.1-29.0%	≧ 29.1%
71-80	≤ 22.5%	22.6-36.5%	36.6-41.5%	≧ 41.6%	≤ 14.0%	14.1-25.0%	25.1-30.0%	≧ 30.1%

### Body Water ( Hydration ) percentage ( % )

Age	Female				Male			
	Low(-)	Normal(0)	Overhigh(+)	High(++)	Low(-)	Normal(0)	Overhigh(+)	High(++)
10-20	≤ 56.0%	56.1-62.0%	62.1-68.0%	≥ 68.1%	≤ 58.0%	58.1-65.0%	65.1-72.0%	≥ 72.1%
21-30	≤ 54.0%	54.1-60.0%	60.1-66.0%	≥ 66.1%	≤ 56.0%	56.1-63.0%	63.1-70.0%	≥ 70.1%
31-40	≤ 51.0%	51.1-57.0%	57.1-63.0%	≥ 63.1%	≤ 54.0%	54.1-61.0%	61.1-68.0%	≥ 68.1%
41-50	≤ 48.0%	48.1-54.0%	54.1-60.0%	≥ 60.1%	≤ 51.0%	51.1-58.0%	58.1-65.0%	≥ 65.1%
51-60	≤ 45.0%	45.1-51.0%	51.1-57.0%	≥ 57.1%	≤ 48.0%	48.1-55.0%	55.1-62.0%	≥ 62.1%
61-70	≤ 42.0%	42.1-48.0%	48.1-54.0%	≥ 54.1%	≤ 45.0%	45.1-52.0%	52.1-59.0%	≥ 59.1%
71-80	≤ 38.0%	38.1-44.0%	44.1-50.0%	≥ 50.1%	≤ 42.0%	42.1-49.0%	49.1-56.0%	≥ 56.1%

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## IMPORTANT INFORMATION - EN

### PACKAGING

In order to respect the environment, packaging material must be disposed of properly in accordance with separate collection.

Check local regulations.

### DISPOSAL

#### DISPOSAL

##### USER INFORMATION

“Implementation of Directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE)”, pertaining to reduced use of hazardous substances in electrical and electronic equipment, as well as to waste disposal.

The symbol of the crossed-out wheeled bin on the equipment or on its packaging indicates that the product must be disposed of separately from other waste at the end of its service life.

The user must therefore take the dismissed equipment to suitable separate collection centres for electrical and electronic waste.

For more details, please contact the appropriate local authority.

Suitable segregated collection of the equipment for subsequent recycling, treatment or environmentally-friendly disposal helps prevent damage to the environment and to human health, and encourages the re-use and/or recycling of the materials that make up the equipment.

Abusive disposal of the product by the user shall result in the application of administrative fines in accordance with the laws in force.

### GUARANTEE

#### Terms and Conditions

**The guarantee is valid for 24 months.**

This guarantee applies only if it has been duly filled in and is submitted with the receipt showing the date of purchase.

The guarantee covers the replacement or repair of parts making up the appliance which were faulty at source due to manufacturing faults.

After the guarantee has expired, the appliance will be repaired against payment.

The manufacturer declines any responsibility for damage to persons, animals or property due to misuse of the appliance and failure to observe the directions contained in the instructions.

#### Limits

All rights under this guarantee and any responsibility on our part will be voided if the appliance has been:

- mishandled by unauthorized persons
- improperly used, stored or transported.

The guarantee does not cover damage to the outward appearance or any other that does not prevent regular operation.

If any faults should be found despite the care taken in selecting the materials and in creating the product, or if any information or advice are required, please contact your local dealer.

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